


21 CFR Part II Compliant

# TOX IVOS<sup>®</sup> II

## Sperm Analyzer

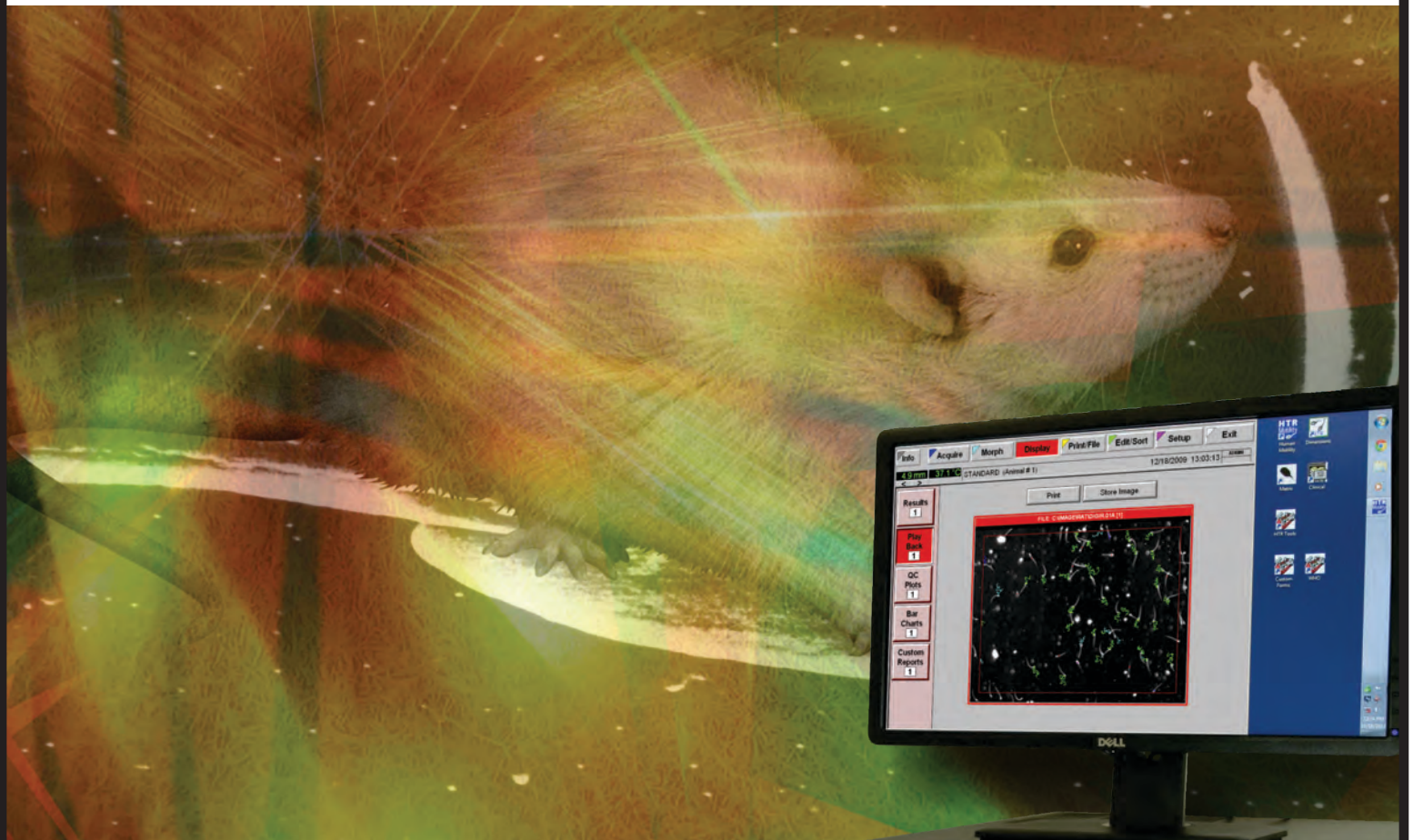

**Fast and Accurate  
Sperm Motility,  
Counts and Viability\***

**Rat, Mouse, Rabbit & More**


\* Mouse and rabbit sperm only. Viability for rat sperm analysis under development.



## Complies with 21 CFR II


### Enhanced Security Features

- ▶ System Administrator assigns User Name, Logon ID, Setup Access, and Temporary Password
- ▶ New user required to create personal password at first logon session
- ▶ Passwords must be 8 - 12 characters in length and contain a mix of letters and numbers
- ▶ User may change password at any time
- ▶ System Administrator may assign a new temporary password for any user at any time
- ▶ Automatic lock out of all users except System Administrator when an invalid Logon ID or password is entered 5 time consecutively

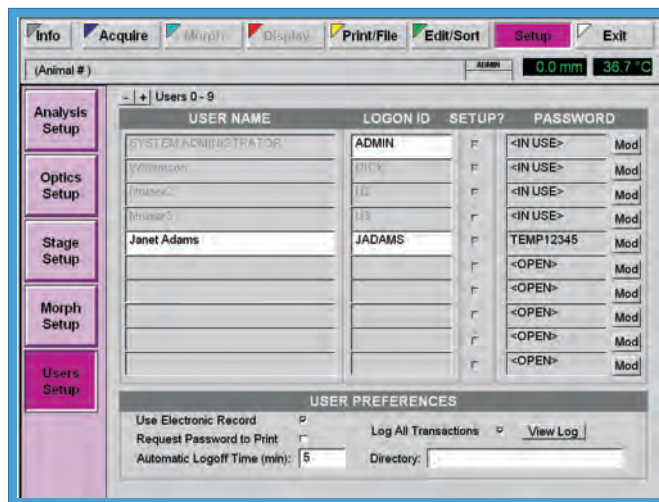

### Electronic Records Keep Data Safe

- ▶ Digital image files stored along with setup files in a single record
- ▶ Previous setup parameters automatically inserted into new Setup 8 to ensure reanalysis with exact settings
- ▶ Any attempts to alter the Electronic Record are detected and reported as corruption when the record is retrieved
- ▶ Printout automatically includes footer containing User Name, User ID, TOX IVOS serial number and software version

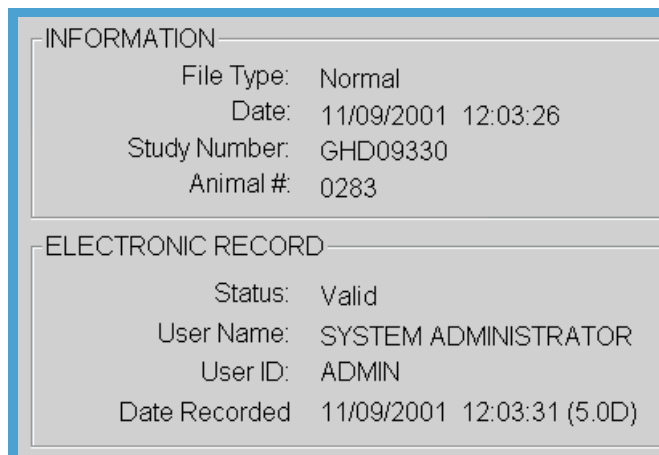

### Log File Tracks Exact Changes

- ▶ Lists User Name, User ID, Date/Time, and Actions performed
- ▶ "Old Setup" and "New Setup" values are shown
- ▶ Multiple attempts at ID or password penetration are logged
- ▶ Log file can be printed by System Administrator

